



Product Catalog

Dental resins, equipments, and accessories

Interactive Catalog

An ISO 13485 Certified Company

[Learn more](#)

TCS is a dental manufacturing company dedicated to providing dental laboratories with exceptional quality materials, equipment, and hands-on training courses that are accessible and reliable. At the core of our brand, we believe that one of the key components to creating successful removable appliances is through forming strong partnerships, because for **TCS**, manufacturing is more than just a business, it is a collectivist culture.

2000

Founded in CA, USA

3

Generations

70+

Countries we sell to

The TCS Experience

Our partnership at a glance

Dynamism

An energetic and passionate team moves forward together and constantly progresses to add value to our customers.

Education

Outstanding technical support and online training designed to grow confidence in creating exceptional TCS dental appliances.

Innovation

A philosophy to continuous improvement guides us to develop better methods, tools, materials, educational programs, etc.

Manufacturing In-House

The majority of our products, including our resins and equipment are manufactured in-house in CA, USA.

Together We Smile - TCS Charitable Events

Our industry is fortunate to be part of a community where giving back literally means putting smiles on people's faces.

Win - Win Philosophy

Your success is our success. We aim for long term relationships based on mutual trust, respect and benefit.

Table of Contents

	Page
Unbreakable TM	6
iFlex TM	10
Karadent TM	14
Millable Discs	18
Injection Units	20
Basic Kits	24
Complete Kits	26
JP90	30
Digital Furnace	32
Waxes / Flasks	33

Table of Contents

	Page
Processing Tools	34
Finishing Tools	35
Polishing Tools	37
TCS Fresh	38
Technical Support	40
Online Training	41
Together We Smile	42
Authorized Distributors	44
Social Media Links	45

UnbreakableTM

Nylon thermoplastic with unbeatable strength



- Extremely resilient, high performance nylon thermoplastic ideal for flexible partials.
- Guaranteed unbreakable due to its unique properties.
- Designed for exceptional memory to hold its form.
- Offers an ideal degree of flexibility.
- A modern fit, lightweight and durable.

[Learn more](#)

Unbreakable™

TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.

Be sure to order the correct diameter cartridge to fit your injection unit.


Millable discs available for CAD/CAM - Refer to pages 18-19.



Shade		25mm	28mm	22mm
Light Pink (Shade 1)	S	4002-01/S1	/	/
	M	4000-01/M1	/	4005-01/M22-1
	L	4005-01/L1	/	4005-01/L22-1
Standard Pink (Shade 2)	S	4002-01/S2	/	/
	M	4000-01/M2	4000-01/MF2	4005-01/M22-2
	L	4005-01/L2	4005-01/LF2	4005-01/L22-2
Light / Dark Pink (Shade 3)	S	4002-01/S3	/	/
	M	4000-01/M3	/	4005-01/M22-3
	L	4005-01/L3	/	4005-01/L22-3
Dark Pink (Shade 4)	S	4002-01/S4	/	/
	M	4000-01/M4	4000-01/MF4	4005-01/M22-4
	L	4005-01/L4	4005-01/LF4	4005-01/L22-4
Natural (Shade N)	S	4002-01/SN	/	/
	M	4000-01/MN	/	4005-01/M22-N
	L	4005-01/LN	/	4005-01/L22-N

Each pack includes the following:
5 delivery bags | **5** patient care instructions | **5** samples of TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:
50 delivery bags and **50** patient care instructions.
*Not applicable to JP Bullets. USA only.**

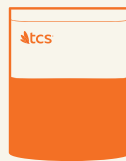

25mm with tab	25 mm with nozzle	1 kilogram	500 grams	JP Bullets
/	/	4003-01/B1	4003-01B/ 1-.5	4002-JP/S1
/	/			/
4005-01/LT25-1	4005-01/LNZ25-1			4005-JP/L1
/	/	4003-01/B2	4003-01B/ 2-.5	4002-JP/S2
/	/			/
4005-01/LT25-2	4005-01/LNZ25-2			4005-JP/L2
/	/	4003-01/B3	4003-01B/ 3-.5	4002-JP/S3
/	/			/
4005-01/LT25-3	4005-01/LNZ25-3			4005-JP/L3
/	/	4003-01/B4	4003-01B/ 4-.5	4002-JP/S4
/	/			/
4005-01/LT25-4	4005-01/LNZ25-4			4005-JP/L4
/	/	4003-01/BN	4003-01B/ N-.5	4002-JP/SN
/	/			/
4005-01/LT25-N	4005-01/LNZ25-N			4005-JP/LN

iFlexTM

Polyolefin thermoplastic with translucent finish



iFlex™

TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.

Be sure to order the correct diameter cartridge to fit your injection unit.

Millable discs available for CAD/CAM - Refer to pages 18-19.



Shade		25mm	28mm	22mm
Light Pink (Shade 1)	S	5002-01/S1	/	/
	M	5000-01/M1	/	5005-01/M22-1
	L	5005-01/L1	/	5005-01/L22-1
Standard Pink (Shade 2)	S	5002-01/S2	/	/
	M	5000-01/M2	5000-01/MF2	5005-01/M22-2
	L	5005-01/L2	5005-01/LF2	5005-01/L22-2
Light / Dark Pink (Shade 3)	S	5002-01/S3	/	/
	M	5000-01/M3	/	5005-01/M22-3
	L	5005-01/L3	/	5005-01/L22-3
Dark Pink (Shade 4)	S	5002-01/S4	/	/
	M	5000-01/M4	5000-01/MF4	5005-01/M22-4
	L	5005-01/L4	5005-01/LF4	5005-01/L22-4
Natural (Shade N)	S	5002-01/SN	/	/
	M	5000-01/MN	/	5005-01/M22-N
	L	5005-01/LN	/	5005-01/L22-N

Each pack includes the following:
5 delivery bags | **5** patient care instructions | **5** samples of TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:
50 delivery bags and **50** patient care instructions.
*Not applicable to JP Bullets. USA only.**



25mm with tab	25 mm with nozzle	1 kilogram	500 grams	JP Bullets
/	/	5003-01/B1	5003-01B/ 1-.5	5005-JP/S1
/	/			/
5005-01/LT25-1	5005-01/LNZ25-1			5005-JP/L1
/	/	5003-01/B2	5003-01B/ 2-.5	5005-JP/S2
/	/			/
5005-01/LT25-2	5005-01/LNZ25-2			5005-JP/L2
/	/	5003-01/B3	5003-01B/ 3-.5	5005-JP/S3
/	/			/
5005-01/LT25-3	5005-01/LNZ25-3			5005-JP/L3
/	/	5003-01/B4	5003-01B/ 4-.5	5005-JP/S4
/	/			/
5005-01/LT25-4	5005-01/LNZ25-4			5005-JP/L4
/	/	5003-01/BN	5003-01B/ N-.5	5005-JP/SN
/	/			/
5005-01/LT25-N	5005-01/LNZ25-N			5005-JP/LN

Karadent™

An upgrade from acrylic without compromising stability



- Semi-rigid, microcrystalline polymer thermoplastic.
- Ideal for full dentures, semi-flexible partial dentures, clear clasps, frameworks, and nightguards.
- An upgrade from traditional acrylic dentures yet as simple to finish and polish.
- Finishes completely transparent in our Clear shade.
- A modern fit, lightweight and durable.

[Learn more](#)

Karadent™





TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.
Be sure to order the correct diameter cartridge to fit your injection unit.



Shade		25mm	28mm	22mm
Standard Pink <i>(Shade 2)</i>	S	2000-01/S2	/	/
	M	2000-01/M2	/	/
	L	2005-01/L2	/	2005-01/L22-2
Dark Pink <i>(Shade 4)</i>	S	2000-01/S4	/	/
	M	2000-01/M4	/	/
	L	2005-01/L4	/	2005-01/L22-4
Natural (clear) <i>(Shade N)</i>	S	2000-01/SN	/	/
	M	2000-01/MN	/	/
	L	2005-01/LN	/	2005-01/L22-N

Each pack includes the following:
5 delivery bags | **5** patient care instructions | **5** samples of TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:
50 delivery bags and **50** patient care instructions.
*Not applicable to JP Bullets. USA only.**


				
25mm with tab	25 mm with nozzle	1 kilogram	500 grams	JP Bullets
/	/	2003-01/B2	2003-01B/ 2-.5	2002-JP/S2
/	/			/
/	/			/
/	/	2003-01/B4	2003-01B/ 4-.5	2002-JP/S4
/	/			/
/	/			/
/	/	2003-01/BN	2003-01B/ N-.5	2002-JP/SN
/	/			/
/	/			/

TCS Millable Discs

Our newest product stands out with **TCS** aesthetics.



[Learn more](#)

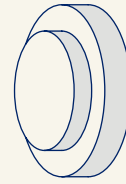

Unbreakable[®]

Shade	98*20mm	98*25mm
Standard Pink (Shade 2)	4000/D98×20-2	4000/D98×25-2
/ Light / Dark Pink (Shade 3)	4000/D98×20-3	4000/D98×25-3
Natural (Shade N)	4000/D98×20-N	4000/D98×25-N

iFlex[®]

Shade	98*20mm	98*25mm
Standard Pink (Shade 2)	5000/D98×20-2	5000/D98×25-2
/ Light / Dark Pink (Shade 3)	5000/D98×20-3	5000/D98×25-3
Natural (Shade N)	5000/D98×20-N	5000/D98×25-N

- Digital technology allows for ultimate precision.
- Appliances come out with a smooth surface making them faster and easier to polish.
- Polishing process is quick (approx. 2 minutes).
- Milling an appliance from a disc is kinder on your tools.

TCS Injection Units

TCS offers a variety of injection units to fit laboratories preferences and budget. View below a few programs available to you.

Equipment Rental Program

Not sure which injection system is right for you?

Get to know TCS injector units before making the investment and avoid the risk of paying restocking fees. Our reasonably priced rental program allows you to test a TCS workflow in your lab for 20 business days.

Get rewarded for smart shopping. Purchase a TCS injection unit (with any authorized distributor) within three months of your rental period and you'll receive free material. USA customers only.

[Click Here to Get A TCS Equipment Rental Quote](#)

SWITCH Promotion Program

Ready to upgrade from another injection brand?

Take advantage of our switch promotions.

[Contact sales@tcsdentalinc.com](mailto:sales@tcsdentalinc.com) for details!


Register your Equipment

Recently purchased a TCS Unit?

Register your equipment to assure your warranty coverage.

[Click Here to Register Your Equipment](#)

TCS Hydraulic Injector

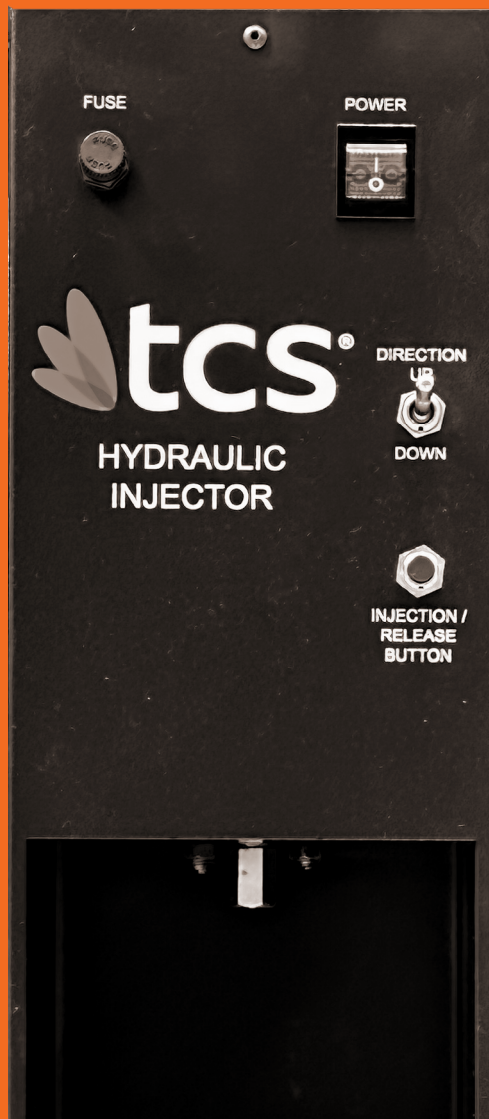

- Equal pressure with every injection
- Minimal maintenance
- No air connections required
- Requires an external heater (TCS Digital Furnace)
- 2 years warranty
- Made in USA

SPECS:

Yields approx. 30 cases per day (*up to 120 cases with additional furnaces*)

H: 22" L: 11" W: 6.5".
50 lbs. (23 kilos)

NOTE:

Used in conjunction with the TCS Digital Furnace. Sold separately.

3504-01	Injector only 110V
3504-02	Injector only 220V

TCS Manual Injector



- Minimal maintenance
- Built to last
- Simple to use
- Requires an external heater (TCS Digital Furnace)
- 2 years warranty
- Made in USA

SPECS:

Yields approx. 30 cases per day (*up to 120 cases with additional furnaces*)

H: 22" L: 22" W: 17"
42 lbs. (19 kilos)


NOTE:

Used in conjunction with the TCS Digital Furnace. Sold separately.

3506-01

Injector only 110V

TCS Automatic Injector

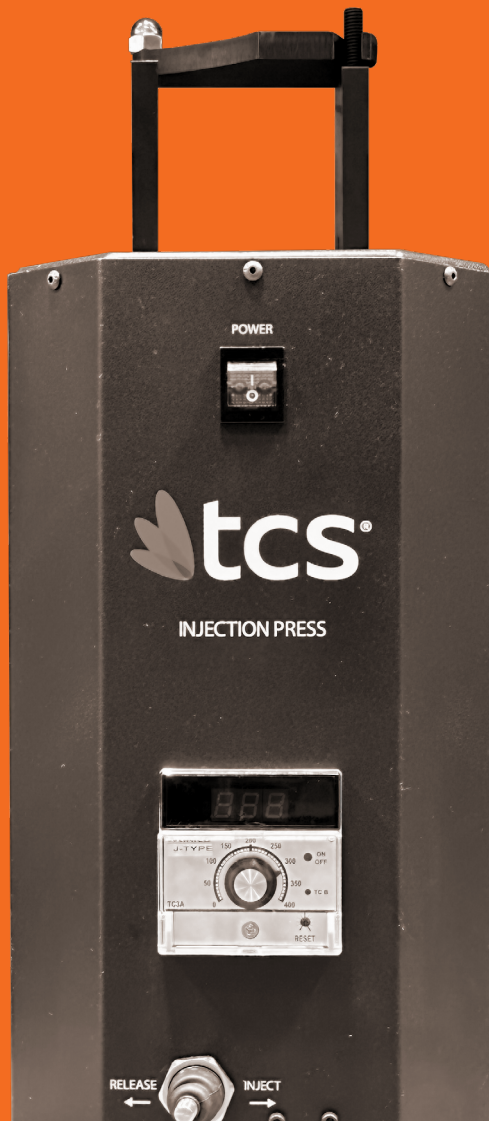

- Compact design with a built-in furnace
- Equal pressure with every injection
- Injects with a touch of a button
- Requires hookup to air pressure (min 140 psi)
- 1 year warranty
- Made in USA

SPECS:

Yields approx. 30 cases per day

H: 23" L: 8" W: 13.5".
40 lbs. (18 kilos).

3525-01	Injector only 110V
3525-02	Injector only 220V

TCS Hydraulic and Manual Injector Basic Kits

Kits include the basic essentials to inject TCS Resins.



1 x Digital Furnace



1 x Cylinder with Handle



1 x Knock Out Rod
1 x Knock Out Stand



1 x Thick Bronze Plug
1 x Thin Bronze Plug



1 x Injection Flask



2 x Cartridges of Each:
Unbreakable™, iFlex™, Karadent™



1 x Mold Release Spray

Hydraulic Basic Kits

3511-01	110V
3511-02	220V

Manual Basic Kits

3508-01	110V
3508-02	220V

TCS Automatic Injector Basic Kits

Kits include the basic essentials to inject TCS Resins.



1 x Injection Flask



2 x Cartridges of Each:
Unbreakable™ and iFlex™



1 x Mold Release Spray

Automatic Air Basic Kits

3528-01	110V
3528-02	220V

TCS Hydraulic and Manual Injector Complete

Includes all the essentials necessary to process TCS Resins.



1 x Digital Furnace



2 x Cylinder with Handle



2 x Injection Flasks



10 pcs. Palate Wax



10 pcs. Lingual Wax



10 pcs. Sprue Wax



3 pcs. 9.5 Sprue Wax



1 x Carving Knife



1 x Cement




1 x EZ Flow Sample



1 x Mold Release Spray



1 x Fusing Liquid

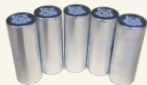

2 pks x Unbreakable™
2 pks x iFlex™
2 pks x Karadent™



1 x Knock Out Rod
1 x Knock Out Stand



2 x Thick Bronze Plug
1 x Thin Bronze Plug



1 x Rubberized Abrasive 304
1 x Rubberized Abrasive 306

Hydraulic Complete Kits

3501-01	110V
3501-02	220V



1 x Hot Air Welder



Twist Drill Samples



TCS Putty (400 gms)



10 x Pink Grinding Wheels



4 x Rag Buffs



1 lb Polycrill



2 x Mini Fuzz Buffs



1 x Turbo Buff



10 x Brown Rubber Wheels



5 x Green Mounted Stones



1 x Pink Zirconia Coarse Bur



2 x Large Rubber Points
2 x Small Rubber Points



Sample Pieces of
Polishing Paste

Manual Complete Kits

3505-01	110V
3505-02	220V

Items in the kits may change without prior notice. We strive for continuous improvement, if we discover a product that works better with our resins, we may switch out an item(s).

TCS Automatic Injector Complete Kits.

Includes all the essentials necessary to process TCS Resins.



2 x Injection Flasks



10 pcs. Palate Wax



10 pcs. Lingual Wax



10 pcs. Sprue Wax



3 pcs. 9.5 Sprue Wax



1 x Carving Knife



1 x Cement



1 x EZ Flow Sample



1 x Mold Release Spray




1 x Fusing Liquid



1 x Hot Air Welder



Twist Drill Samples

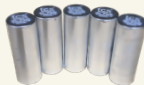

3 pks x Unbreakable™
3 pks x iFlex™



2 x Large Rubber Points
2 x Small Rubber Points



1 x Rubberized Abrasive 304
1 x Rubberized Abrasive 306

Automatic Complete Kits

3526-01	110V
3526-02	220V



TCS Putty (400 gms)



10 x Pink Grinding Wheels



10 x Brown Rubber Wheels



2 x Mini Fuzz Buffs



1 x Pink Zirconia Coarse Bur



5 x Green Mounted Stones



4 x Rag Buffs



1 x Turbo Buff



1 lb Polycril



Sample Pieces of
Polishing Paste

Note: It is not recommended to inject Karadent™ by tcs® using the Automatic Air Injector.

Items in the kits may change without prior notice. We strive for continuous improvement, if we discover a product that works better with our resins, we may switch out an item(s).

JP90 Handheld Injector



- 1-3 tooth unilaterals
- Up to 5 consecutive tooth partials/saddles, saddles for combination cases
- Repairs
- 2 years warranty
- Made in USA

SPECS:

H: 6" L: 2-1/4" W: 11"
2 lbs. (1 kilo)

7000-JP

JP90 Injector only



JP Flask

Fabricated from high precision stainless steel to withstand warpage under extreme heat.

7005-JP



Metal Sprues

Used to create a gateway for injection.

7014-JP



Matrix Pusher

Used to facilitate putty application.

7010-JP



TCS Putty

Used to create a matrix when injecting with the JP90 Handheld Injector. Can also be used to invest for clean injections.

4296-01

Base & Catalyst
400gms ea.

JP90 Starter Kits



1 x JP90 Handheld
Injector



JP Bullets (Small)
1 pk x Unbreakable™
1 pk x iFlex™



JP Bullets (Large)
1 pk x Unbreakable™
1 pk x iFlex™



JP flask



TCS Putty (400 gms)



2 x Waxed Cases for
Practice



1 x Matrix Pusher
5 x Metal Sprues

JP90 Starter Kits

7001-JPU	JP Starter Kit
7020-JPU	JP Starter Kit with Digital Furnace

Digital Furnace & Accessories



- Simple and practical
- Melts thermoplastic resin
- Compact design
- Reliable / low maintenance
- 2 years warranty
- Made in USA

SPECS:
H: 6" L: 1.5" W: 11"
1.3 lbs. (0.6 kilos)

3512-01	110V
3512-02	220V



Melting Cylinder with Handle

Fabricated from high precision stainless steel to withstand warpage under extreme heat.

3703-01



Knock Out Stand and Rod

Steel tools used for purging the melting cylinder from the remaining cartridge and bronze plug after injection.

3705-01	Knock Out Stand
3707-01	Knock Out Rod



Bronze Plugs

Thick bronze plug is used when melting medium and small cartridges in the furnace. Thin bronze plug is used when melting large cartridges in the furnace.

3709-01	Bronze Plug Thin
3712-01	Bronze Plug Thick

Waxes / Flasks



Palate Wax

Save time waxing with pre-formed palates created with the recommended thickness of a TCS palate.

5015-01 **80/box**



Lingual Wax

Pre-formed lingual bars are ideal for clasps, flange, and lingual connector areas.

5020-01 **1 lb**



Sprue Wax

Easy to work with, does not get brittle. Ideal diameter for spruing. 5 gauge.

5025-01 **1 lb**



Sprue Wax 9.5

Connect your sprue wax to the 9.5 Sprue Wax only at the tip of the flask opening to avoid unwanted plaster sticking to it.

5027-01	5 pcs
5026-01	1 lb (~30 pcs)
Each piece yields ~10 injections.	



Injection Flask

Milled from one piece of aluminum block. Extra smooth finish facilitates the de-flasking process. Includes a special knock-out button to simplify case removal and prolong the life of the flask.

3701-01



Pour Flask

Used for acrylic pour technique and/or to duplicate models. Manufactured in light-weight aluminum. Easy to assemble, disassemble and clean.

3650-01



TCS Carving Knife

Excellent for removing extended fibers that sometimes appear when processing TCS partials. Also used to define the necklines of the TCS restorations when necessary.

4201-01	
4202-01	Replacement Blades (12 pcs)

Processing Tools



TCS Cement

To secure teeth in the stone prior to injection. Mild enough to keep teeth from bonding to the stone during injection.

3802-01 **0.5 oz**



Fusing Liquid

For repairs and additions to existing TCS restorations. The Fusing Liquid is ONLY used when flasking and injecting are necessary.

3812-01 **0.5 oz**



Mold Release Spray

Colorless and non-staining multi-purpose lubricating spray, used to ease the purging of the furnace accessories and TCS cartridges during injection process.

3814-01 **12 oz**



Hot Air Welder

Hand torch in a slim design for minor repairs or additions to TCS restorations. Refillable with butane. (Butane not included).

- * 6 month warranty
- * proof of purchase (invoice) is required

3714-01



Twist Drills

To create the diatoric retention in acrylic teeth necessary for fabricating TCS restorations.

005	5005-01	6 pcs
006	5006-01	6 pcs
009	5009-01	6 pcs
012	5012-01	6 pcs



Empty Cartridges

Empty cartridges are sold in packs of 50

	25mm	28mm	22mm
S	3003-01S	/	/
M	3003-01M	3003-01MF	3003-01M22
L	3003-01L	3003-01LF	3003-01L22



EZ Flow

Separating medium for processing thermoplastics. Eliminates stone sticking to the partial, leaving it shiny and clean.

4217-01	16 oz
4218-01	1 gallon

Finishing Tools



Pink Grinding Wheels

Excellent and economical product for the finishing stage.

4205-01 100/box



Brown Rubber Wheels

Use prior to pumicing in order to obtain a high shine. It is crucial that the entire partial is rubberized and smoothened.

4503-01 100/box



Green Mounted Stones

Used during the finishing stage. Cuts thermoplastics moderately, and leaves a smooth surface. Ideal for minor chairside adjustments.

4501-01 100/box



Knock Down Wheel Chuck

Mandrel used with the White Knock Down Wheel and Rubberized Abrasives.

4212-01



Rubberized Abrasives

Long-lasting wheel for smoothening the partial during the finishing stage.

4304-01	Thin
4306-01	Thick



Brown Rubber Points

Ideal to eliminate scratches and smoothen the surface.

S	4513-01	12 pcs
L	4514-01	12 pcs



Mini Fuzz Buff

Use after the Brown Rubber Wheels to further smoothen partial prior to pumicing.

4222-01 12 pcs



Carbide Bur

Specifically designed to cut thermoplastics. Long-lasting with high cutting efficiency.

50251G-01

Finishing Tools / Polishing Tools



Pink Zirconia Fine Bur

Use during the finishing stage to leave the surface smooth and prepped for the Brown Rubber Wheel. Also, ideal for chairside adjustments. Max RPM= 12,000. Use with a hand piece.

5035Z



Pink Zirconia Coarse Bur

Use during the finishing stage to cut excess material and trim appliance to the desired length. Gross reduction. Max RPM= 16,000. Use with a hand piece

5045Z



TCS Embroidered Apron

Adjustable apron with pockets and embroidered TCS logo.

Apron01

1 pcs



Polycril

An advanced formula for pumicing acrylics and thermoplastics at high-speed.

4298-01

1 lb



Turbo Buff

Pumicing buff with silicone cotton cloth. Use for pumicing only.

4208-01



Rag Buff

Stitched only in the center for the extra softness needed for polishing TCS restorations.

4207-01



Mini Rag Buff

Recommended to dentists for polishing touch-ups after minor adjustments on TCS appliances. Use with Green Hi-Shine polishing paste.

4221-01

12 pcs

Polishing Tools



Brown Tripoli

Use after pumicing to prepare appliance for polishing. Effective for all TCS resins. Use with a dry Rag Buff.


4206-01 2 lbs



iShine

Use after Brown Tripoli to polish surface of iFlex™ partials (and other polyolefin thermoplastics). Use with a dry Rag Buff.

4220-01 10.75 oz

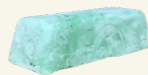

Green Hi-Shine

Use after Brown Tripoli to polish surface of Karadent™ and Unbreakable™ partials. Use with a dry Rag Buff.

4209-02 1 lb



Dr. Adjustment Kit.*

Ideal for quick and easy chairside adjustments.

4014-02



**Sample Model
Unbreakable™**

An ideal marketing tool for your laboratory.

4011-01U	Upper Partial
4011-01L	Lower Partial



**Sample Model:
Karadent™**

An ideal marketing tool for your laboratory.

4016-01U	Upper Full Denture
4017-01L	Lower Full Denture
4018-01U	Upper Clear Framework
4019-01U	Upper Partial Denture
4019-01L	Lower Partial Framework



**Sample Model:
iFlex™**

An ideal marketing tool for your laboratory.

4009-01U	Upper Partial
4009-01L	Lower Partial





A clean appliance is vital to patients' health and the appearance and longevity of dental appliances. TCS Concentrated Dental Appliance Cleanser is a safe solution for quick and effective results.



TCS Concentrated Dental Appliance Cleaner

Professional strength cleanser. Deep cleans, removes stains and odors. Ideal for most types of removable dental appliances.

4505-02

24 pcs
6 month supply

TCS Dental Appliance Cleaner

TCS Dental Appliance Cleaner sample packets. Box of 75 packets (1-week supply ea.).

4505FR-SMP

75 pcs

TCS Dental Appliance Brush

Resilient soft bristle brush to thoroughly clean dental appliances without scratching the surface.

4506FR-01

TCS includes patient care instructions and a sample pack of denture cleaner with each cartridge.*
Include them with the finished case so that the patient can maintain a clean and vibrant appliance.

**USA only*

TCS Support

We value your success!

Quality tools and support help improve your laboratory's productivity and quality when fabricating a TCS Dental Appliance.

Text Support

Register your phone for quick text support:

Text your name, name of laboratory, city, and state to (562) 241-1225

Request Promo Material

Grow your removable department. Request digital web-ads and printed promotional material. *Ship to US only.

Order Now

Online Training

Sign up for an eLearning experience. TCS provides thorough follow-along training for producing TCS Flexible Partial The TCS Way. 2 CDT credits will be awarded. More information may be found on the following page.

Start Today





Online Training

[Start Today](#)

Covers All the Basics

Covers design, block out, setup, waxing, processing, finishing and polishing.

Basic Training Kit

Kits are available for USA customers only and includes all the basic tools needed to follow along with the course.

Test your Knowledge

Each chapter includes a quiz so you can test how much you've learned.

Certificate of Completion

Receive a digital certificate of completion upon completing the course.

Earn CDT Credits

The training course gets you 2 Scientific CDT credits!

Pricing

\$80

Complete: Includes a Welcome Box containing a variety of tools to follow along with each lesson. Ships to USA only.

\$42

Basic: Access to online training only. Use coupon code BASIC at checkout.

Free

Basic training with purchase of any TCS Injector Kit. Contact sales@tcsdentalinc.com for free access within 30 days of your purchase.

Together We Smile

The Together We Smile-TCS Gives Back Charitable Event is an outreach program that travels to remote areas across the globe to deliver partials and dentures to people in need. Although oral care is available to everyone, those who live on low incomes have no real access to proper oral care and treatment. TCS launched Together We Smile to serve those communities and bring happy smiles whenever we can!



TCS Together We Smile was launched in 2018 with its first mission in the Dominican Republic. Partnering with Jack Meibaum and his charity, Mission Para Los Niños, TCS set up a pop-up dental lab to treat and fabricate partials and dentures for over 60 patients. This initiative was repeated in 2019, furthering our commitment to restoring smiles.

In 2022, our team traveled to La Guajira, Colombia, followed by San Andrés, Colombia, in 2023, and Thailand in 2024. Each mission is dedicated to providing as many patients as possible with dentures and partials within a one to one-and-a-half-week timeframe.

Most recently, in early January 2025, our team journeyed to Sangolquí, Ecuador, where we had the privilege of transforming the lives of over 46 patients through our dental care services, creating over 80 appliances in just four days. The work is demanding, with long hours and tireless effort, but the joy and gratitude from our patients make it all worthwhile.

Carol Barros, from Ecuador, along with her family and a group of local volunteers, carried out a patient pre-selection process several months in advance. This preparation ensured that when the full TWS team arrived, the workflow would be more efficient. Mr. Oscar Barros' laboratory (Carol's father) was made available to us throughout the four-day event.

The impression-taking station was set up in an unoccupied community family home. Portable dental chairs were rented, and the space was fully adapted to accommodate all patients. The location was ideal due to its accessibility for those receiving treatment.

The Together We Smile outreach program thrives today thanks to the collaboration of local dentists, technicians, generous dental suppliers, and our incredible volunteers. By working together, we continue to provide life-changing dental solutions to individuals in need, restoring confidence and happiness—one smile at a time.



[Read more](#)

TCS Authorized Distributors



Alvy Dental Supply
(800) 426-3626



American Dent-All, Inc.
(818) 662-0618




Benco Dental
(800) 462-3626



Dentorium
(800) 847-4209



Exceed Dental
(833) 539-2333

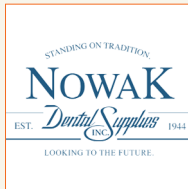

**Nowak
Dental Supplies**
(800) 654-7623




**Pan America
Dental Supply**
(305) 266-3630



**Pearson
Dental Supply**
(800) 336-8256




RMH3 Dental Supply
(571) 269-2776

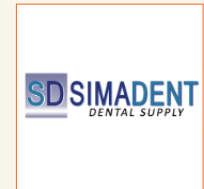

Simadent USA
(305) 827-2811



Unicore Dental Supply
(800) 581-9356

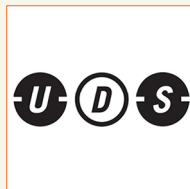

Uhler Dental Supply
(800) 937-3753



Zahn Dental
(800) 496-9500



TCS Dental Inc.



@tcsdentalinc



TCS Dental Inc.



TCS Dental Inc.

2619 Lime Ave.,
Signal Hill, CA 90755, USA

(866) 426-2970 | (562) 426-2970
www.tcsdentalinc.com